

Quality inside

Neo-Safe® T CU 380

Intrauterinpešsar zur Empfängnisverhütung
Intrauterine Contraceptive Device
Dispositif Contraceptif Intra-Utérin
Dispositivo contraccettivo intrauterino

1 I.U.D.

CU-Safe® T 300

Intrauterinpešsar zur Empfängnisverhütung
Intrauterine Contraceptive Device
Dispositif Contraceptif Intra-Utérin
Dispositivo contraccettivo intrauterino

1 I.U.D.

Multi-Safe® CU 375

Intrauterinpešsar zur Empfängnisverhütung
Intrauterine Contraceptive Device
Dispositif Contraceptif Intra-Utérin
Dispositivo contraccettivo intrauterino


1 I.U.D.

T-Safe® CU 380A QL

Intrauterinpešsar zur Empfängnisverhütung
Intrauterine Contraceptive Device
Dispositif Contraceptif Intra-Utérin
Dispositivo contraccettivo intrauterino

1 I.U.D.

EURIM
PHARM

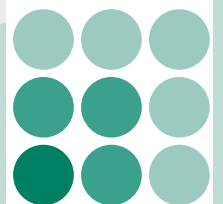

IUDs according to international standards

All IUDs, distributed by EurimPharm, are generics of the classic copper IUDs that have been on the market for many years.

The body of the IUD is made of Polyethylene with copper wire wound around the stem. The addition of Barium Sulphate enables radiological proofing.

The total copper surface area varies between 300 mm² and 380 mm² allowing a different in-situ-time of the IUDs (between 5 and 10 years).

The copper has a purity level of 99,99 %. Undesirable allergic reactions can therefore be excluded.

The return thread at the bottom end of the body is composed of tear-proof Nylon.

EurimPharm Intrauterine Devices –

T-Safe[®]

Optimal, secure fit in the uterus: minimal expulsion rate.

Available in two versions:

T-Safe CU 380A

- Copper surface: 380 mm²
- In-situ-time: 10 years
- IUD with the lowest Pearl Index

T-Safe CU 380A QL


- Identical characteristics to T-Safe CU 380A
- Easy and quick loading of the IUD due to the new technology

CU-Safe[®]

Same design as „Flexi T 300“
Easy, painfree insertion and removal
Secure double nylon thread system

CU-Safe T 300

- Copper surface: 300 mm²
- In-situ-time: 5 years
- Optimally suited for nullipara
- Also suited for uni- and multipara with smaller uteri
- Stem length: 29,2 mm

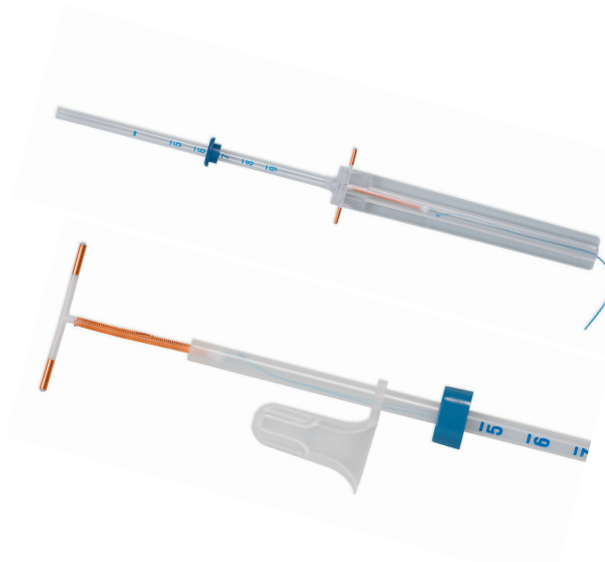

Quality advantage through EurimPharm IUDs

All products carry the internationally valid CE-sign and are certified in accordance with the thorough Quality Management System (actually valid ISO 13485) that is subject to numerous controls through the national surveillance authorities.

Complete physician's instruction and patient information according to new standard (ISO 7439:2015 Copper-bearing contraceptive intrauterine devices — Requirements and tests)

Loading capsule with the T-Safe CU 380A

→ Easy insertion of the side arms into the insertion tube

Quickload insertion system with the T-Safe CU 380A QL

Tear-proof and yet soft, pliant nylon threads and safe, firm sitting copper cover

→ Minimal irritation in the vagina

Secure, tear-proof sterile packaging (Tyvek)

High quality at low prices

Multi-Safe®

Form and composition same as the long-year tried and trusted best seller „Multiload“.
One-handed insertion technique.

Available in two versions:

Multi-Safe CU 375

- Copper surface: 375 mm²
- In-situ-time: 5 years

Multi-Safe CU 375 short

- Identical characteristics to Multi-Safe CU 375
- Shorter stem length (29,4 mm). Ideally suited for a shorter uterus length

Neo-Safe®

Same design as „Nova T“
Easy, painfree insertion and removal


Available in two versions:

Neo-Safe T CU 380

- Copper surface: 380 mm²
- In-situ-time: 5 years
- Low Pearl Index
- Low PID (Pelvic inflammatory disease)

Neo-Safe T CU 380 Mini

- Identical characteristics to Neo-Safe T CU 380
- Optimally suited for nulipara

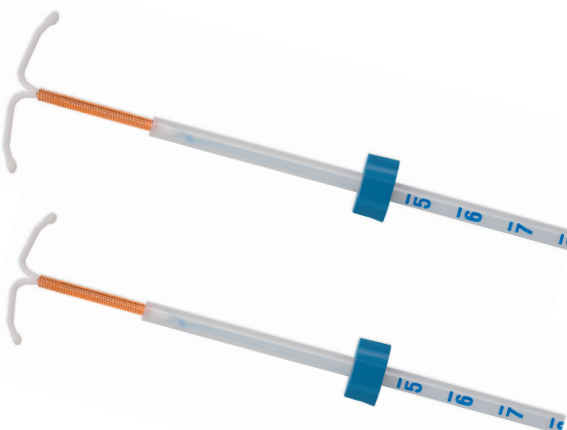

EurimPharm IUDs

High quality at low prices

Today EurimPharm offers all current IUD models, each with various copper surface areas. The selection of the most suitable IUD depends in the main on individual patient wishes and requirements.

Useful Information about EurimPharm IUDs

| EurimPharm IUD | Width (mm) | Length (mm) | Maximum Insertion Time in the Body | Recommended Uterus Length |
|-------------------------|------------|-------------|------------------------------------|---------------------------|
| Multi-Safe CU 375 | 19,5 | 35,2 | 5 years | 6 – 9 cm |
| Multi-Safe CU 375 short | 19,5 | 29,4 | 5 years | 5 – 8 cm |
| T-Safe CU 380A | 31,9 | 35,9 | 10 years | 6 – 9 cm |
| T-Safe CU 380A QL | 31,9 | 35,9 | 10 years | 6 – 9 cm |
| Neo-Safe T CU 380 | 31,8 | 31,9 | 5 years | 6 – 9 cm |
| Neo-Safe T CU 380 Mini | 24,0 | 30,0 | 5 years | 5 – 8 cm |
| CU-Safe T 300 | 22,9 | 29,2 | 5 years | 5 – 8 cm |

Whichever EurimPharm IUD you select, you are choosing
high quality at low prices

Magnetic Resonanz Therapie (MRT/MRI)

MR-examinations in a strength of up to 3.0 Tesla in women who have a copper intrauterine device for contraception (IUD) in place, does not influence neither the results of the MR-examination nor the functional efficacy of the IUD.

For more information please see www.eurim-iud.de.

CE
0120

Contact:

EurimPharm Vertriebs-GmbH & Co KG
Michael-Walz-Gasse 18d
A-5020 Salzburg

Ms. Elisabeth Adomaitis
Mr. Thomas Valerius
Phone: +49 8654 7707-327 or -146
Telefax: +49 8654 7707-189
Email: elisabeth.adomaitis@eurim.de
thomas.valerius@eurim.de
Internet: www.eurim-iud.de

